



### Micro current beauty instrument

### Instruction manual

#### Attention

1. Pregnant women, minors under 18, do not use
2. Please do not use for patients with heart disease or diabetes, or for people with pacemaker, internal defibrillator, nerve stimulator, insulin pump and other electronic instruments;
3. Do not use metal implants (such as gold wire) within the scope of use;
4. Do not use if you have cancer, suspected cancer or precancerous lesions, especially skin precancerous nevus (carnivorous nevus);
5. Serious concurrent diseases such as heart disease, epilepsy, AIDS and immune system damage do not use;
6. Do not use if there are skin diseases such as sores, psoriasis, eczema, rash, recurrent herpes simplex, etc;
7. Please do not use if you have had a scar, abnormal wound healing or fragile skin;
8. Laser skin rejuvenation, chemical peeling and wound healing, please avoid or use after complete healing, if there are more medical problems, please consult a doctor;
9. Within the scope of use, temporary fillers (hyaluronic acid, water light needle, protein thread) and permanent fillers (prosthesis nose augmentation and chin augmentation) have been implanted for more than three months and can be used after consultation with your doctor;
10. Do not use for people with metal allergy;
11. Do not use for people with sensitive skin;
12. In the process of using, please stop using immediately in case of any discomfort.

#### Maintenance and service

1. Protect the instrument from moisture. When cleaning, prevent water from entering USB port and charging port. Do not put the instrument into dishwasher or microwave oven.
2. Please turn off the instrument before cleaning.
3. Do not use sharp or corrosive, volatile diluent cleaner.
4. You can use a soft cloth or make-up cotton that is wet, supercooled or lukewarm.
5. Do not wash the instrument directly with water, because it will cause water damage to the instrument.

1

#### Product Brief

Congratulations on receiving the micro current beauty instrument  
To make sure you get results immediately, read the entire user manual.  
This beauty instrument is a small portable micro current facial care device, which is used for facial feeling  
Situation improvement.  
The instrument will gently stimulate the large surface area of the face and provide 5 minutes of face lifting service,  
Improve your appearance. The device transmits low-order microcurrent pulses to the chin, cheekbones and forehead  
The upper part.  
Double sphere design can be used to make the best contact with the face of various shapes and sizes. Allow users  
Adjust the output level of micro current to achieve personalized comfort.

Not more than twice a day.

2

#### Matters needing attention

- ▲ Supervise children to make sure they don't play with the instrument.
- ▲ Put the packaging materials such as plastic bags in places that children can't touch.
- ▲ Do not use the instrument near the bathtub, washbasin or water environment.
- ▲ Please use a charger that meets the requirements of this instrument to charge.
- ▲ If the instrument is damaged, in order to avoid danger, it must be repaired by the manufacturer and professional maintenance personnel.
- ▲ Do not drop or knock the instrument. If the instrument is damaged due to dropping or collision, do not use it.
- ▲ If the instrument accidentally falls into the water, please shut down immediately and do not try to start it, and contact the factory after-sales service.
- ▲ Please follow the operation instructions and do not use the machine for other purposes.
- ▲ When the instrument is abnormally hot, please stop using and contact the factory for after-sales service.
- ▲ Because the instrument is equipped with polymer battery, do not place the instrument in a high temperature environment or in the sun.
- ▲ Do not connect the USB port and the charging port of the instrument together with the charging line to avoid accidents.


#### After-sale service

1. This instrument is only suitable for personal care, not for professional use.
2. We provide a one-year warranty for our products, which starts from the date of purchase.
3. During the warranty period, we will replace parts or the whole product as the case may be due to defects in materials or production process.
4. The warranty claim must be submitted to the authorized repair station or factory for after-sales service together with the certificate (purchase invoice) that can prove the purchase date.

3

#### Component description

1. Power on / off and strength control up / down button
2. Output level indicator
3. Micro current treatment head
- 4.USB interface
5. Speaker (not shown):
  - A. It beeps when the power is on.
  - B. It beeps when the power is off.
  - C. When the intensity is adjusted, it will beep.
  - D. Double beeps when moving to a new position on the skin.
  - E. When the device is low on power, two series of beeps indicate that the battery needs to be charged.

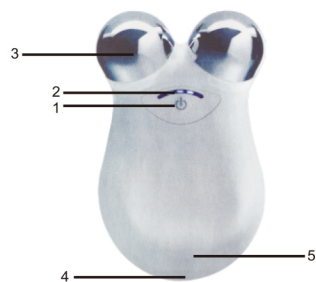

#### Using steps

1. Please follow all the steps in the "before using the device" section.
2. Please clean your face with oil-free detergent. Men need to shave before care because facial hair may interfere with the conductivity of the sphere.
3. please apply gelatin or emulsion to clean, dry skin that needs care. Do this after making sure there is no allergic reaction to gel or emulsion.

4

5. The company does not guarantee
  - Damage caused by improper use, such as using the wrong power supply, inserting the product into an unsuitable socket, artificial damage, etc.
  - Normal use causes product aging and wear, but has little effect on the use function of the product.
  - Product problems caused by improper operation or non professional maintenance or use of non company accessories for other purposes. (after sales service does not include damage caused by the above reasons)

#### Technical parameters and packaging

Product Name: micro current beauty instrument  
Working voltage: 2.4V  
Output voltage: 5V  
Product size: 106mm \* 65mm \* 38mm  
Power: 0.15w  
Net weight: 128G  
Battery: 2.4v/800mah polymer battery  
Packing list: one host, one USB charging line, one cloth bag, one manual

Note: the product pictures and the introduction of relevant product features as well as the explanation of some terms shown in this manual can only be used as a reference. When users use the product, the material object shall prevail. In view of the consistent policy of continuous improvement of the product, the company reserves the right to change the product specifications without prior notice.

5

4. Please gradually press the on / off button to increase or decrease the micro current output and adjust the intensity level. When you switch between three intensity levels, the instrument will beep and the number of blue lights will change.  
Strength grade:  
Minimum intensity: left 2 lights on  
Medium intensity: the middle three lights are on  
Maximum intensity: two lights on the right
5. Please contact the instrument directly with the area to be cared, see "care area" below. Keep both spheres on your face at the same time. Note: micro current beauty instrument will not contract muscles. When used correctly, the current is gentle, very soothing, and gentle.
6. please do not dry gel or emulsion during the nursing process. The gel or emulsion can be reacted by coating more gel or emulsion in the nursing area.
7. after finishing nursing, wipe the gel or lotion with damp cloth.

Nursing area <https://www.youtube.com/shorts/7RAspZHF9I8>



6